

Incubators

HettCube 600 | 600 R

Efficient use of valuable space

Building on the base of the 400 model, the HettCube 600 has a maximum validated usable space of 351 L. This model is capable of holding up to 1,800 petri dishes.

This unit is available with refrigeration and a temperature range from 0 °C to +65 °C (HettCube 600 R).



Highlights

- Maximal validated usable space on a small footprint
- One-hand operation without additional inner door
- 4.3 inch touch display for intuitive operation
- Automatic door closure with magnetic seals, door orientation can be easily changed
- Low noise level of ≤ 44 dB(A)
- TÜV certified
- Factory certificate HettCert with 9 measuring points according to DIN 12880:2007-5 standard
- HettCube without cooling requires only 1 °K (1 °C) above stable ambient temperature
- No interior fan
- Model variants with IVDR approval including DIN 62304 standard
- Inclusive: 3 standard shelves and 1 HTS shelf with telescopic rail

Additionally for cooled models:

- Temperature drop with holiday function
- Temperature selection monitor class 3.2
- Easy integration of external devices via the switch-board (optional)

Feature

- At a glance display of: settings, future programs and target/actual values, temperature drops, events
 - real-time calendar
 - More individual settings: Flexible door alarm, asymmetric tolerance band
 - Three different power failure scenarios
 - Logbook: records of every door opening, events, alarms and operating hours
 - Language options EN, DE, FR, ES
 - Wide range of program functions: Start after time, start after temperature, standstill segment
 - Programs and temperature settings transferable to other HettCubes without additional software
 - Interior made of high-quality stainless steel
- Special accessories:
- Shelves and racks for storage of samples
 - Shelf for Petri Dishes
 - Tightly welded HTS drawers (various heights)

Technical data

	HettCube 600 incubators	HettCube 600 R cooled incubators
Temperature range	1 K above ambient temperature up to +65 °C	0 °C up to +65 °C
Exterior dimensions (without bushing and door handle) W x D x H in mm	710 x 825 x 1990	710 x 825 x 1990
Interior dimensions W x D x H in mm	535 x 690 x 1415	535 x 690 x 1415
Internal volume in liters	520	520
Validated usable volume in liters	351	351
Percentage validated usable volume / internal volume	67 %	67 %
Footprint in m ²	0,6	0,6
Weight in kg	164	175
Number of trays provided as standard	4 (3 standard + 1 HTS)	4 (3 standard + 1 HTS)
Temperature fluctuation at +37 °C	± 0.1 K	± 0.1 K
Temperature uniformity at +37 °C	± 0.2 K	± 0.2 K
Temperature uniformity at +25 °C	± 0.1 K	± 0.1 K
Recovery time after door has been opened for 30 s at +37 °C	≤ 5,5 min	≤ 5,5 min
Energy consumption at +37 °C	0.046 kWh/h	0.046 kWh/h
Noise level	≤ 41 dB (A)	≤ 44 dB (A)
Power supply	220 – 240 V 1 ~ / 50 – 60 Hz	220 – 240 V 1 ~ / 50 – 60 Hz
Cat. No.	66000	66005
Non-IvD version	66001	66006
Other voltages		
100-120 V 1 ~ / 50 – 60 Hz	66000-01	66005-01



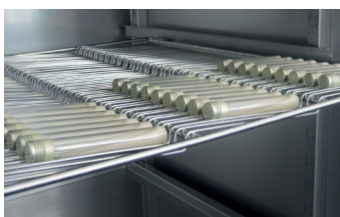
Watch our animation video
of the HettCube on
YouTube.



options and accessories



Switchboard

Rack for Loewenstein
application

Rack for petri dishes

	Cat. No.
Shelf (Set) Made of stainless steel, with standard rails, max. load (kg): 50	60001
Shelf (HTS-Set) Made of stainless steel, with telescopic rails, extendable up to 70 %, max. load (kg): 40	60031
Drawer (HTS-Set), High 30 mm Made of stainless steel, with telescopic rails, extendable up to 70 %, tightly welded, max. load (kg): 40	60024
Drawer (HTS-Set), High 65 mm Made of stainless steel, with telescopic rails, extendable up to 70 %, tightly welded, max. load (kg): 40	60025
Drawer (HTS-Set), High 105 mm Made of stainless steel, with telescopic rails, extendable up to 70 %, tightly welded, max. load (kg): 40	60026
Rack (HTS-Set) For Petri dishes, stainless steel, with telescopic rails, extendable up to 70 %, Petri dishes Ø (mm): 90, max. load (pcs): 90	60038
Rack (Set) For Petri dishes, stainless steel, with standard rails, Petri dishes Ø (mm): 90, max. load (pcs): 90	60039
Rack For Petri dishes, stainless steel, Petri dishes Ø (mm): 90, max. load (pcs): 90	60040
Rack (HTS-Set) For inclined storage of cultures (Loewenstein), stainless steel, with telescopic rails, extendable up to 70 %, inclination angle: 5°, tube Ø (mm): 15-20, max. load (pcs): 81 tube	60036
Rack (Set) For inclined storage of cultures (Loewenstein), stainless steel, with standard rails, inclination angle: 5°, tube Ø (mm): 15-20, max. load (pcs): 81 tubes	60037
Rack For inclined storage of cultures (Loewenstein), stainless steel, inclination angle: 5°, tube Ø (mm): 15-20, max. load (pcs): 81 tubes	60041
Frame L, 16-place Made of stainless steel, for inclined storage of cultures, tube Ø (mm): 15-20, Tube length (mm): 100-125, inclination angle 5° or 20°	60027
Frame XL, 16-place Made of stainless steel, for inclined storage of cultures, tube Ø (mm): 15-20, Tube length (mm): 126-170, inclination angle 5° or 20°	60028
Switchboard 4-fold socket strip, as a unit controllable via display, on the back of the device	60521
Independent PT 100 sensor For independent temperature measurement, four-wire system, temperature values output with analogue output 4-20 mA on the back of the device	60503
Passive dehumidification For individual or timed opening of a dehumidification module via the touchscreen	60042
Service: Assembly of the Stacking kit for the HettCube 200 200 R	60043
Service: Changing the door hinge	60044
Glass door All-glass outer door, for HettCube 200 200 R	60030
Glass door All-glass outer door, for HettCube 400 400 R	60029
Glass door All-glass outer door, for HettCube 600 600 R	60013
Access port Ø (mm): 22, Foam stoppers	60006
Access port Ø (mm): 42, Foam stoppers	60007
Access port Ø (mm): 67, Foam stoppers	60008
Stacking kit For safe stacking of two HettCubes 200 200 R	60009
Rolling cabinet Lockable, with one drawer, incl. lockable castors, WxDxH (mm): 770x800x550, for HettCube 200 200 R	60010
USB Port Lock (Set) For securing the USB-A interface. Set consisting of 10 securing clips and 1 USB key tool.	60525
Parameter for continuous cooling operation For storing samples below 15 °C for more than 2 weeks	60526

— extensive standard equipment

4.3 inch touchscreen:

- Target / actual display
- Setting accuracy 0.1 ° C
- Real-time calendar
- Timer
- Language options (English, German, French and Spanish)
- Temperature diagram in 3 zoom levels (up to 4 weeks)
- Power failure scenarios
- Door alarm individually adjustable
- Log (door openings, alarms and operating hours)
- PIN lock
- Up to 99 program functions (Start/Stop, Period, Timer at start, Timer at temperature, ...)
- Temperature selection monitor class 3.1 for all models

Additionally with refrigerated devices

- Temperature reduction with holiday function
- Temperature selection monitor class 3.2 for all refrigerated models
- Programming external devices via switchboard (optional)

Control panel and door lockable simultaneously

USB service interface

Bushing on rear panel Ø 42 mm

Potential-free alarm output

Interior of high-quality stainless steel (W-St 1.4301 (ASTM 304))

In addition you get for free ...

3 standard shelves of stainless steel

1 HTS shelf of stainless steel with telescopic rails

Factory certificate (HettCert) – 9 points measurement analog to DIN 12880: 2007-05

Compensation for unneeded standard shelves*

* Choose other shelves / drawers as standard, we credit you the price of unused standard inserts.

— Certifications / Registrations

Product certification:



Product registration:



Company certifications:



— Dimensions

