




No rights can be derived from these images. These images refer to a specific model. This model is high (not the standard height).

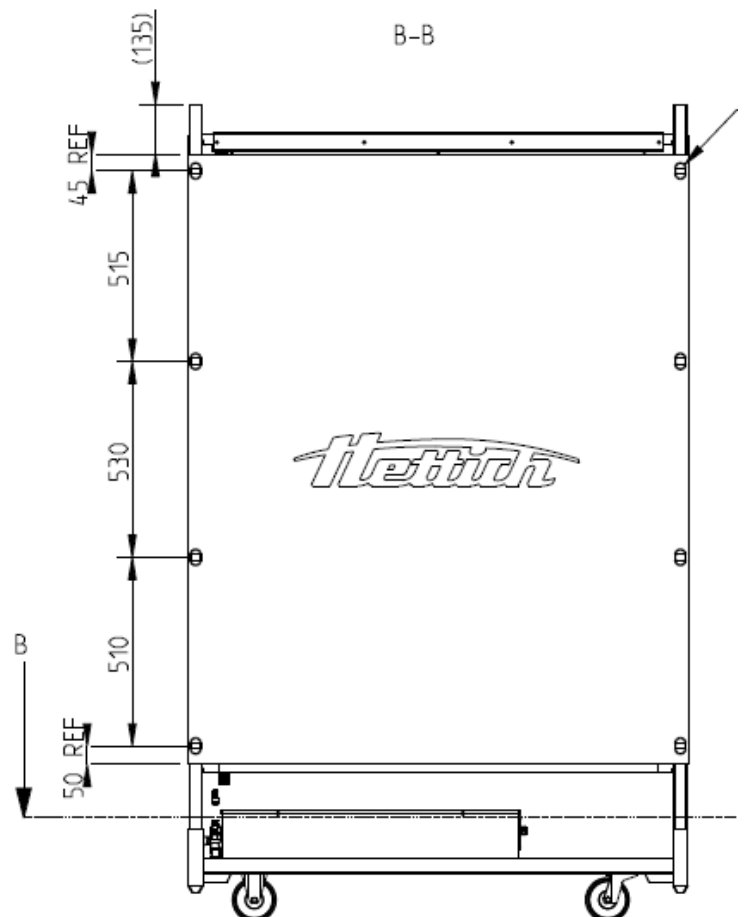

TROLLEY / DANISH CART

Making your day to day work so much more easier and flexible. The Danish cart on its own or as part of your Hettich Benelux climate room.

The Hettich Benelux, Danish carts are particularly popular in horticulture. Thanks to their ideal dimensions, manoeuvrability and adjustable platforms, these carts are not only ideal for transport. They are also easy to use as mobile racks in, for example, a climate chamber. The cart can be equipped with the 600 or 1100 µmol LEDs and an optional watering system. To protect your plants from damaging, we also offer screens.

Specifications

Standard Model

External dimensions d x w x h (mm)	1350 x 565 x 1900 mm
Internal dimensions d x w x h (mm)	1270 x 545 x 1735mm
Material	Coated galvanised steel
Gross weight	15,5 kg
Wheels	2x swivel, 2x fixed castors
Number of trays	max. 3
LED (dimmable)	600 or 1100 µmol Incl. FarRed (dimmable)
Type of light	constant current
Homogeneity light 600 µmol	at 15 cm
Homogeneity light 1100 µmol	at 50 cm
Max. weight per shelf	100 kg
Dimmable lighting	10 -100 %
Colour	White
Optional	
Watering system	per tray, with quick valve
Screens	Breathable, grey, protection cloths
Extra height model dimensions	1350 x 565 x 2100 mm
Technical requirements	
Electrical supply	230 V, 50/60 Hz

Type

Type	Art. No.
Standard model, 600 µmol	8740.0600
Standard model, 1100 µmol	8740.1100

Custom-built Controlled Environment Room HSP50

In our factory, we build custom environment rooms for pharmaceutical stability or plant growth. This is an example. This room was dismantled for delivery to the customer and later reassembled on site. This custom-built room is designed to accommodate 16 Danish trolleys, each with its own independent light control and watering system. Also equipped with a touch screen for easy management.

With dimensions of 6.77 m x 4.88 m x 3.10 m and internal volume of 50 m³, it offers ample space for your needs.

The possibilities for a custom product of Hettich Benelux are endless.

